

Code No. 27417

Human Thioredoxin Assay Kit - IBL

INTRODUCTION

Human thioredoxin is a 12 kDa* protein which consists of 105 amino acids and regulates the redox status in living bodies by reactions of two vicinal cysteines in its active site.

It is known that thioredoxins are induced by stress such as ultraviolet, radiation, oxidant and viral infection, and act as a biological defense factor against stress, showing antioxidant effects and regulating activities of intracellular signaling molecules and transcription factors.

Therefore, measuring thioredoxin is expected to monitor body's stress status, and the study to improve the stress-resistance has been proceeding.

PRINCIPLE

This kit is a solid phase sandwich ELISA using 2 kinds of high specific antibodies. Tetra Methyl Benzidine (TMB) is used as a coloring agent (Chromogen). The strength of coloring is proportional to the quantities of Human Thioredoxin1.

MEASUREMENT RANGE

3.91 ~ 250 ng/mL

INTENDED USE

- This IBL's assay kit is capable for the quantitative determination of Human Thioredoxin1 in serum, EDTA plasma, urine and cell culture supernatant.
- In the case of measuring culture medium with FCS, we recommend to set the negative medium control as a test sample blank.

KIT COMPONENT

1	Precoated plate :		
	Anti-Human Thioredoxin (5t1) Mouse IgG MoAb Affinity Purify	96Well x 1	
2	Labeled antibody Conc. :		
	(30X) HRP conjugated Anti- Human Thioredoxin (N) Rabbit IgG Fab' Affinity Purify	0.4mL x 1	
3	Standard :	Recombinant Human Thioredoxin	1.0mL x 2
4	EIA buffer		30mL x 1
5	Solution for Labeled antibody		12mL x 1
6	Chromogen : TMB solution		15mL x 1
7	Stop solution		12mL x 1
8	Wash buffer Conc.		50mL x 1

OPERATION MANUAL

1. Materials needed but not supplied

- Plate reader (450nm)
- Graduated cylinder and beaker
- Incubator (37°C ± 1°C)
- Paper towel
- Washing bottle for precoated plate
- Disposable test tube for "2, Labeled antibody Conc." and "6, Chromogen"
- Micropipette and tip
- Deionized water
- Graph paper (log/log)
- Tube for dilution of Standard


2. Preparation

- 1) Preparation of wash buffer
 "8, Wash buffer Conc." is a concentrated (40X) buffer. Adjust the temperature of "8, Washing buffer Conc." to room temperature and then, mix it gently and completely before use. Dilute 50 mL of "8, Wash buffer Conc." with 1,950 mL of deionized water and mix it. This is the wash buffer for use. This prepared wash buffer shall be stored in refrigerator and used within 2 weeks after dilution.
- 2) Preparation of Labeled antibody
 "2, Labeled antibody Conc." is a concentrated (30X). Dilute "2, Labeled antibody Conc." with "5, Solution for Labeled antibody" in 30 times according to required quantity into a disposable test tube. Use this resulting solution as Labeled antibody.
 Example)
 In case you use one strip (8 well), the required quantity of Labeled antibody is 800 µL. (Dilute 30 µL of "2, Labeled antibody Conc." with 870 µL of "5, Solution for Labeled antibody" and mix it. And use the resulting solution by 100 µL in each well.)
This operation should be done just before the application of Labeled antibody.
 The remaining "2, Labeled antibody Conc." should be stored at 4°C in firmly sealed vial.
- 3) Preparation of Standard
 Put just 1.0 mL of deionized water into the vial of "3, Standard" and mix it gently and completely. This solution is 500 ng/mL Human Thioredoxin standard.
- 4) Dilution of Standard
 Prepare 8 tubes for dilution of "3, Standard". Put 230 µL each of "4, EIA buffer" into the tube.
 Specify the following concentration of each tube."

Tube-1	250 ng/mL
Tube-2	125 ng/mL
Tube-3	62.5 ng/mL
Tube-4	31.25 ng/mL
Tube-5	15.63 ng/mL
Tube-6	7.81 ng/mL
Tube-7	3.91 ng/mL
Tube-8	0 ng/mL (Test Sample Blank)

Put 230 µL of Standard solution into tube-1 and mix it gently. Then, put 230 µL of tube-1 mixture into tube-2. Dilute two times standard solution in series to set up 7 points of diluted standard between 250 ng/mL and 3.91 ng/mL. Tube-8 is the test sample blank as 0 ng/mL.

See following picture.

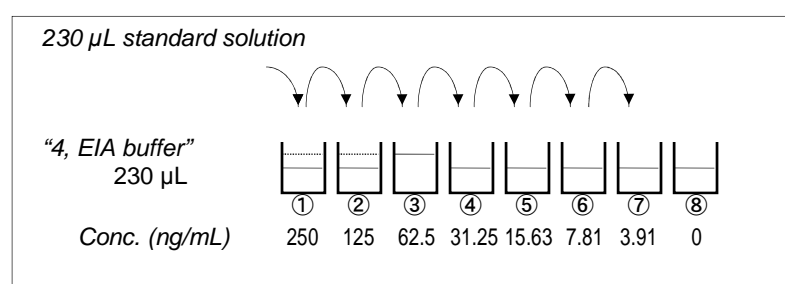

5) Dilution of test sample

Test sample should be diluted with "4, EIA buffer" as necessary.

If the concentration of Human Thioredoxin in samples may not be estimated in advance, the pre-assay with several different dilutions will be recommended to determine the proper dilution of samples.

3. Measurement procedure

All reagents shall be brought to room temperature approximately 30 minutes before use. Then mix it gently and completely before use. Make sure of no change in quality of the reagents. Standard curve shall be prepared simultaneously with the measurement of test samples.

Reagents	Test Sample	Standard	Test Sample Blank	Reagent Blank
	Test sample 100 µL	Diluted standard (Tube 1-7) 100 µL	EIA buffer (Tube-8) 100 µL	EIA buffer 100 µL
Incubation for 60 minutes at 37°C with plate lid				
4 times (wash buffer more than 350 µL)				
Labeled Antibody	100 µL	100 µL	100 µL	-
Incubation for 30 minutes at 37°C with plate lid				
5 times (wash buffer more than 350 µL)				
Chromogen	100 µL	100 µL	100 µL	100 µL
Incubation for 30 minutes at room temperature (shielded)				
Stop solution	100 µL	100 µL	100 µL	100 µL
Read the plate at 450nm against a Reagent Blank within 30 minutes after addition of Stop solution.				

- 1) Determine wells for reagent blank. Put 100 µL each of "4, EIA buffer" into the wells.
- 2) Determine wells for test sample blank, test sample and diluted standard. Then, put 100 µL each of test sample blank (tube-8), test sample and dilutions of standard (tube-1-7) into the appropriate wells.
- 3) Incubate the precoated plate for 60 minutes at 37°C after covering it with plate lid.
- 4) Wash the plate with the prepared wash buffer and remove all liquid.
- 5) Pipette 100 µL of labeled antibody solution into the wells of test samples, diluted standard and test sample blank.
- 6) Incubate the precoated plate for 30 minutes at 37°C after covering it with plate lid.
- 7) Wash the plate with the prepared wash buffer and remove all liquid.
- 8) Take the required quantity of "6, Chromogen" into a disposable test tube. Then, pipette 100 µL from the test tube into the wells. Please do not return the rest of the test tube to "6, Chromogen" bottle to avoid contamination.
- 9) Incubate the precoated plate for 30 minutes at room temperature in the dark. The liquid will turn blue by addition of "6, Chromogen".
- 10) Pipette 100 µL of "7, Stop solution" into the wells. Mix the liquid by tapping the side of precoated plate. The liquid will turn yellow by addition of "7, Stop solution".
- 11) Remove any dirt or drop of water on the bottom of the precoated plate and confirm there is no bubble on the surface of the liquid. Then, run the plate reader and conduct measurement at 450 nm against a reagent blank. The measurement shall be done within 30 minutes after addition of "7, Stop solution".

SPECIAL ATTENTION


- 1) Test samples should be measured soon after collection. For the storage of test samples, store them frozen and do not repeat freeze/thaw cycles. Thaw the test samples at a low temperature and mix them completely before measurement.
- 2) Test samples should be diluted with "4, EIA buffer", if the need arises.
- 3) Duplicate measurement of test samples and standard is recommended.
- 4) Use test samples in neutral pH range. The contaminations of organic solvent may affect the measurement.
- 5) Use only wash buffer contained in this kit for washing the precoated plate. Insufficient washing may lead to the failure in measurement.
- 6) Remove the wash buffer completely by tapping the precoated plate on paper towel. Do not wipe wells with paper towel.
- 7) "6, Chromogen" should be stored in the dark due to its sensitivity against light. "6, Chromogen" should be avoided contact with metals.
- 8) Measurement should be done within 30 minutes after addition of "7, Stop solution".

CALCULATION OF TEST RESULT

Subtract the absorbance of test sample blank from all data, including standards and unknown samples before plotting. Plot the subtracted absorbance of the standards against the standard concentration on log-log graph paper. Draw the best smooth curve through these points to construct the standard curve. Read the concentration for unknown samples from the standard curve.

Example of standard curve

Conc. (ng/mL)	Absorbance (450nm)
250	1.595
125	1.236
62.5	0.868
31.25	0.593
15.63	0.407
7.81	0.279
3.91	0.199
0 (Test Sample Blank)	0.091

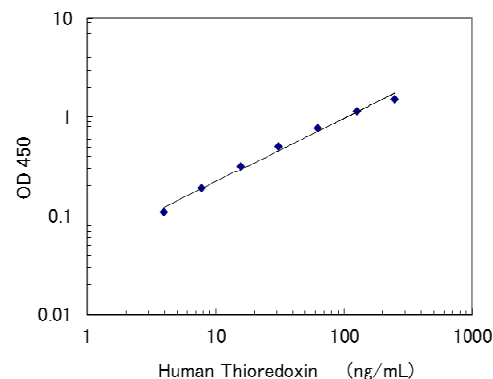

* The typical standard curve is shown above. This curve can not be used to derive test results. Please run a standard curve for each assay.

PERFORMANCE CHARACTERISTICS

1. Titer Assay (Samples with standard added are used.)

Specimen	Titer (X)	Measurement Value (ng/mL)	Theoretical Value (ng/mL)	%
10%FCS added RPMI-1640	2	160.11	155.19	103.2
	4	81.76	75.35	108.5
	8	40.72	35.81	113.7
Human Serum	2	107.06	125.00	85.6
	4	54.26	62.50	86.8
	8	29.51	31.25	94.4
Human Plasma (EDTA)	2	107.60	125.00	86.1
	4	54.76	62.50	87.6
	8	30.29	31.25	96.9
Human Urine	2	127.86	125.00	102.3
	4	70.61	62.50	113.0
	8	35.11	31.25	112.4

2. Added Recovery Assay

Specimen	Theoretical Value (ng/mL)	Measurement Value (ng/mL)	%
10%FCS added RPMI-1640 (x2)	157.37	153.30	97.4
	94.87	108.05	113.9
	63.62	73.57	115.6
Without FCS RPMI-1640	125.00	115.90	92.7
	62.50	53.59	85.7
	31.25	25.66	82.1
Human Serum (x4)	33.25	28.14	84.6
	17.63	14.95	84.8
	9.81	8.09	82.5
Human Plasma (EDTA) (x4)	31.25	25.86	82.8
	15.63	12.74	81.5
	7.81	6.37	81.6
Human Urine (x2)	32.02	27.20	84.9
	16.40	14.92	91.0
	8.58	7.87	91.7

3. Intra - Assay

Measurement Value (ng/mL)	SD value	CV value (%)	n
110.44	10.00	9.1	20
30.75	2.88	9.4	20
10.65	0.96	9.0	20

4. Inter - Assay

Measurement Value (ng/mL)	SD value	CV value (%)	n
109.08	7.92	7.3	6
34.26	4.48	13.1	6
11.71	1.25	10.7	6

5. Specificity

Compound	Cross Reactivity
Human Thioredoxin1	100 %
Human Thioredoxin2	≤0.1 %
Human VEGF	≤0.1 %
Human VEGF-C	≤0.1 %
Human Osteopontin	≤0.1 %
Human Leptin	≤0.1 %

6. Sensitivity

0.43 ng/mL

The sensitivity for this kit was determined using the guidelines under the National Committee for Clinical Laboratory Standards (NCCLS) Evaluation Protocols. (National Committee for Clinical Laboratory Standards Evaluation Protocols, SC1, (1989) Villanova, PA: NCCLS.)

PRECAUTION FOR INTENDED USE AND/OR HANDLING

- All reagents should be stored at 2 - 8°C. All reagents shall be brought to room temperature approximately 30 minutes before use.
- "3, Standard" is lyophilized products. Be careful to open this vial.
- "7, Stop solution" is a strong acid substance. Therefore, be careful not to have your skin and clothes contact "7, Stop solution" and pay attention to the disposal of "7, Stop solution".
- Dispose used materials after rinsing them with large quantity of water.
- Precipitation may occur in "2, Labeled antibody Conc.", however, there is no problem in the performance.
- Wash hands after handling reagents.
- Do not mix the reagents with the reagents from a different lot or kit.
- Do not use expired reagents.
- This kit is for research purpose only. Do not use for clinical diagnosis.

STORAGE AND THE TERM OF VALIDITY

Storage Condition : 2 - 8°C
The expiry date is specified on outer box.

Version 4.

November 2016 *

Made in Japan.