

Code No. 27410

Mouse angiopoietin-like 3 Assay Kit - IBL

INTRODUCTION

Angiopoietin-related protein 3 (ANGPTL3) is secreted by the liver and is a hormonelike protein involved in lipid metabolism and angiogenesis. ANGPTL3 is structurally similar to the angiopoietin family, which is involved in angiogenesis. It possesses a coiled-coil domain and fibrinogen-like domain, but no binding to the Tie family has been observed. ANGPTL3 has been found to be under gene expression regulation by LXR, and it has attracted interest as a new secretion factor.

PRINCIPLE

This kit is a solid phase sandwich ELISA using 2 kinds of high specific antibodies. Tetra Methyl Benzidine (TMB) is used as coloring agent (Chromogen). The strength of coloring is in proportion to the quantities of Mouse angiopoietin-like 3.

MEASUREMENT RANGE

0.31 \sim 20 ng/mL

INTENDED USE

The IBL's Mouse angiopoietin-like 3 Assay Kit is a complete kit for the quantitative determination of Mouse angiopoietin-like 3 in serum, EDTA-plasma and cell culture media.

KIT COMPONENT

1	Precoated plate	: Anti-Mouse angiopoietin-like 3 Rabbit IgG Affinity Purify	96Well x 1
2	Labeled antibody	Conc.	
	: (30X) HRP conjugate	ed Anti-Mouse angiopoietin-like 3 Rabbit IgG Fab' Affinity Purify	0.4mL x 1
3	Standard	: Recombinant Mouse angiopoietin-like 3	0.5mL x 2
4	EIA buffer		30mL x 1
5	Solution for Label	ed antibody*	12mL x 1
6	Chromogen	: TMB solution	15mL x 1
7	Stop solution*		12mL x 1
8	Wash buffer Conc	· ·	50mL x 1

OPERATION MANUAL

1. Materials needed but not supplied

- Micropipette and tip
- Plate reader (450nm)Graduated cylinder and beaker
 - Deionized water
 Graph paper (log/log)
- Refrigerator (as 4°C)
- Tube for dilution of Standard
- Paper towel
- Incubator (37℃±1℃)
- · Washing bottle for precoated plate
- Disposable test tube for "2, Labeled antibody Conc." and "6, Chromogen"

2. Preparation

1) Preparation of wash buffer

"8, Wash buffer Conc." is a concentrated (X40) buffer. The temperature of "8, Wash buffer Conc." shall be adjusted to room temperature and then, mix it gently and completely before use. Dilute 50mL of "8, Wash buffer Conc." with 1,950mL of deionized water and mix it. This is the wash buffer for use. This prepared wash buffer shall be stored in refrigerator and used within 2 weeks after dilution.

2) Preparation of Labeled antibody

"2, Labeled antibody Conc." is a concentrated (X30). Dilute "2, Labeled antibody Conc." with "5, Solution for Labeled antibody" in 30 times according to required quantity into a disposable test tube. Use this resulting solution as Labeled antibody.

Example)


In case you use one slit (8 well), the required quantity of Labeled antibody is 800 μ L. (Dilute 30 μ L of "2, Labeled antibody Conc." with 870 μ L of "5, Solution for Labeled antibody" and mix it. And use the resulting solution by 100 μ L in each well.)

This operation should be done just before the application of Labeled antibody. The remaining "2, Labeled antibody Conc." should be stored at 4°C in firmly sealed vial.

- 3) Preparation of Standard
- Put just 0.5 mL of deionized water into the vial of "3, Standard" and mix it gently and completely. This solution is 40 ng/mL Mouse angiopoietin like 3 standard. 4) Dilution of Standard
- Prepare 8 tubes for dilution of "3, Standard". Put 230 μ L each of "4, EIA buffer" into the tube.

On a sife the falles do not a set of a set of a set to be

See following picture.

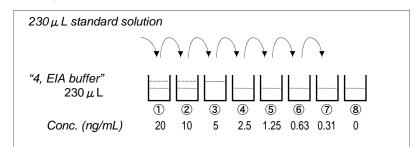

5) Dilution of test sample

Test sample should be diluted with "4, EIA buffer" as the need arises. If the concentration of Mouse angiopoietin-like 3 in samples may not be estimated in advance, the pre-assay with several different dilutions will be recommended to determine the proper dilution of samples.

3. Measurement procedure

All reagents shall be brought to room temperature approximately 30 minutes before use. Then mix it gently and completely before use. Confirm no change in quality of the reagents. Standard curve shall be prepared simultaneously with the measurement of test samples.

	Test Sample	Standard	Test Sample Blank	Reagent Blank	
Reagents	Test sample 100 μ L	Diluted standard (Tube 1~7) 100 µ L	EIA buffer (Tube-8) 100 μ L	EIA buffer 100 μ L	
	Incubation for 60 minutes at 37°C with plate lid				
	4 times (wash buffer more than 350 μL)*				
Labeled Antibody	100 µ L	100 µ L	100 µ L	-	
Incubation for 30 minutes at 4°C with plate lid					
	5 times (wash buffer more than 350 µL)*				
Chromogen	100 μ L	100 μ L	100 μ L	100 μ L	
Incubation for 30 minutes at room temperature (shielded)					
Stop solution	100 μ L	100 μ L	100 µ L	100 μ L	
Read the plate at 450nm against a Reagent Blank within 30 minutes after application of Stop solution.					

- 1) Determine wells for reagent blank. Put $100 \,\mu$ L each of "4, EIA buffer" into the wells.
- 2) Determine wells for test sample blank, test sample and diluted standard. Then, put 100 μ L each of test sample blank (tube-8), test sample and dilutions of standard (tube-1~7) into the appropriate wells.
- 3) Incubate the precoated plate for 60 minutes at 37°C after covering it with plate lid.
- 4) Wash the plate with the prepared wash buffer and remove all liquid.*
- 5) Pipette $100 \,\mu$ L of labeled antibody solution into the wells of test samples, diluted standard and test sample blank.
- Incubate the precoated plate for 30 minutes at 4°C after covering it with plate lid.
- 7) Wash the plate with the prepared wash buffer and remove all liquid.*
- 8) "6, Chromogen" should be taken the required quantity into a disposable test tube. Then, pipette $100 \,\mu$ L from the test tube into the wells. Please avoid to return the rest of test tube into "6, Chromogen" bottle due to avoid to cause of contamination.
- 9) Incubate the precoated plate for 30 minutes at room temperature in the dark. The liquid will turn blue by the addition of "6, Chromogen".
- 10) Pipette $100 \,\mu$ L of "7, Stop solution" into the wells. Mix the liquid by tapping the side of precoated plate. The liquid will turn yellow by the addition of "7, Stop solution".
- Remove any dirt or drop of water on the bottom of the precoated plate and confirm there is no bubble on the surface of the liquid. Then, run the plate reader and conduct measurement at 450nm. The measurement shall be done within 30minutes after the addition of "7, Stop solution".

SPECIAL ATTENTION

1) Test samples should be measured soon after the collection. In case of the storage of test samples, they should be stored under frozen conditions and do

Specify the following concentration of each tube.

(Test Sample Blank)
(

Put $230 \,\mu$ L of Standard solution into tube-1 and mix it gently. Then, put $230 \,\mu$ L of tube-1 mixture into tube-2. Dilute two times standard solution in series to set up 7 points of diluted standard between 20 ng/mL and 0.31 ng/mL. Tube-8 is the test sample blank as 0 ng/mL.

not repeat freeze/thaw cycles. Thaw the test samples at low temperature and mix them completely before measurement.

- 2) Test samples should be diluted with "4, EIA buffer", if the need arises.
- 3) The measurement of test samples and standard in duplicate is recommended.
- 4) Use test samples in neutral pH range. The contaminations of organic solvent may affect the measurement.
- 5) Use only wash buffer contained in this kit for washing the precoated plate. Insufficient washing may lead to the failure in measurement.
- 6) Remove the wash buffer completely by tapping the precoated plate on paper towel.

Do not wipe wells with paper towel.

- 7) "6, Chromogen" should be stored in the dark due to its sensitivity against light."6, Chromogen" should be avoided contact with metals.
- 8) Measurement should be done within 30 minutes after addition of "7, Stop solution".



CALCULATION OF TEST RESULT PERFORMANCE CHARACTERISTICS

1. Titer Assay (Samples with standard added are used.)

Specimen	Titer (X)	Measurement Value (ng/mL)	Theoretical Value (ng/mL)	%
10%FCS	4	2.26	2.50	90.4
added	8	1.22	1.25	97.6
RPMI-1640	16	0.66	0.63	104.8
	2	13.17	14.22	92.6
Mouse Serum (Balb/c)	4	8.08	7.21	112.1
(Balb/0)	8	4.58	3.62	126.5
Mouse Plasma	2	11.97	11.29	106.0
(EDTA)	4	6.76	6.11	110.6
(Balb/c)	8	3.67	3.16	116.1

2. Added Recovery Assay

Specimen	Theoretical Value (ng/mL)	Measurement Value (ng/mL)	%
	5.00	5.14	102.8
10%FCS added RPMI-1640 (x4)	2.50	2.50	100.0
	1.25	1.24	99.2
	12.27	13.08	106.6
Mouse Serum (Balb/c) (x8)	11.65	12.15	104.3
	11.33	11.44	101.0
Mouse Plasma	11.39	11.84	104.0
(EDTA) (Balb/c) (x8)	8.89	9.53	107.2
	7.64	8.40	109.9

3. Intra - Assay

Measurement Value (ng/mL)	SD value	CV value (%)	n
7.75	0.28	3.6	24
2.04	0.13	6.4	24
0.61	0.04	6.6	24

4. Inter - Assay

Measurement Value (ng/mL)	SD value	CV value (%)	n
7.54	0.53	7.0	34
2.01	0.13	6.5	34
0.61	0.04	6.6	34

6. Sensitivity

0.04 ng/mL

The sensitivity for this kit was determined using the guidelines under the National Committee for Clinical Laboratory Standards (NCCLS) Evaluation Protocols. (National Committee for Clinical Laboratory Standards Evaluation Protocols, SC1, (1989) Villanova, PA: NCCLS.

PRECAUTION FOR INTENDED USE AND/OR HANDLING

- All reagents should be stored at 2~8°C. All reagents shall be brought to room 1. temperature approximately 30 minutes before use.
- 2.
- "3, Standard" is lyophilized products. Be careful to open this vial."7, Stop solution" is a strong acid substance. Therefore, be careful not to contact 3. your skin and clothes with "7, Stop solution" and pay attention to the disposal of "7, Stop solution".
- 4. Dispose used materials after rinsing them with large quantity of water.
- The precipitation may grow in "2, Labeled antibody Conc.", however, there is no 5. problem in the performance.
 - Wash hands after handling reagents. 6.
 - Do not mix the reagents with the reagents from different lot or different kit. 7.
 - 8. Do not use the reagents expired.
 - 9. This kit is for research purpose only. Do not use for clinical diagnosis.

STORAGE AND THE TERM OF VALIDITY

Storage Condition : 2 - 8 °C The expiry date is specified on outer box.

REFERENCE

- 1. Koishi R, Ando Y, Ono M, Shimamura M, Yasumo H, Fujiwara T, Horikoshi H, Furukawa H. Angptl3 regulates lipid metabolism in mice. Nat Genet. 2002 Feb; 30(2): 151-7.
- Camenisch G, Pisabarro MT, Sherman D, Kowalski J, Nagel M, Hass P, Xie MH, 2. Gurney A, Bodary S, Liang XH, Clark K, Beresini M, Ferrara N, Gerber HP. ANGPTL3 stimulates endothelial cell adhesion and migration via integrin alpha v beta 3 and induces blood vessel formation in vivo. J Biol Chem. 2002 May 10; 277(19): 17281-90.

Version 2.

September 2016*

Immuno-Biological Laboratories Co., Ltd.