

Code No. 27758

Human Napsin A Assay Kit - IBL

INTRODUCTION

Napsin A (NAPSA) is a pepsin-like aspartic proteinase that belongs to the peptidase A1 family and plays a role in pneumocyte surfactant processing. Napsin A is expressed as a single chain protein with the molecular weight of approximately 38 kDa. Immunohistochemical studies revealed high expression levels of Napsin A in human lung and kidney but low expression in spleen. Napsin A is expressed in type II pneumocytes and in adenocarcinomas of lung. The high specificity expression of napsin A in adenocarcinomas of lung is useful to distinguish primary lung adenocarcinomas from adenocarcinomas of other organs. This ELISA kit can measure human Napsin A.

PRINCIPLE

This kit is a solid phase sandwich ELISA using 2 kinds of highly specific antibodies. Tetra Methyl Benzidine (TMB) is used as a coloring agent (Chromogen). The strength of coloring is proportional to the quantities of Human Napsin A.

MEASUREMENT RANGE

3.91 - 250 ng/mL

INTENDED USE

This IBL's assay kit is capable for the quantitative determination human Napsin A in serum, urine and cell culture supernatant.

KIT COMPONENT

- Precoated plate 1
- : Anti-Human Napsin A (61A2) Mouse IgG MoAb Affinity Purify 96Well x 1 2 Labeled antibody Conc. : (30X) HRP conjugated Anti- Human Napsin A (TMU-Ad02) Mouse IgG MoAb Fab' Affinity Purify 0.4mL x 1 3 Standard : Reconbinant Human Napsin A 0.5mL x 2 EIA buffer 30mL x 1 4 12mL x 1 Solution for Labeled antibody* 5 15mL x 1 6 Chromogen : TMB solution 12mL x 1 Stop solution* 7 Wash buffer Conc.* 8 50mL x 1

OPERATION MANUAL

1. Materials needed but not supplied

- Plate reader (450nm)
 - Graduated cylinder and beaker · Deionized water
- Refrigerator (as 4°C)
- · Graph paper (log/log) Tube for dilution of Standard

· Micropipette and tip

- · Paper towel · Washing bottle for precoated plate
- · Disposable test tube for "2, Labeled antibody Conc." and "6, Chromogen"

2. Preparation

- 1) Preparation of wash buffer
 - "8, Wash buffer Conc." is a concentrated (40X) buffer. Adjust the temperature of "8, Washing buffer Conc." to room temperature and then, mix it gently and completely before use. Dilute 50 mL of "8, Wash buffer Conc." with 1,950 mL of deionized water and mix it. This is the wash buffer for use. This prepared wash buffer shall be stored in refrigerator and used within 2 weeks after dilution.
- 2) Preparation of Labeled antibody

"2, Labeled antibody Conc." is a concentrated (30X). Dilute "2, Labeled antibody Conc." with "5, Solution for Labeled antibody" in 30 times according to required quantity into a disposable test tube. Use this resulting solution as Labeled antibody.

Example)

In case you use one strip (8 well), the required quantity of Labeled antibody is 800 µL. (Dilute 30 µL of "2, Labeled antibody Conc." with 870 µL of "5, Solution for Labeled antibody" and mix it. And use the resulting solution by 100 μ L in each well.)


This operation should be done just before the application of Labeled antibody. The remaining "2, Labeled antibody Conc." should be stored at 4°C in firmly sealed vial.

Preparation of Standard 3)

Put just 0.5 mL of deionized water into the vial of "3, Standard" and mix it gently and completely. This solution is 500 ng/mL Human Napsin A standard. **Dilution of Standard**

Prepare 8 tubes for dilution of "3, Standard". Put 230 µL each of "4, EIA buffer" into the tube.

See following picture.

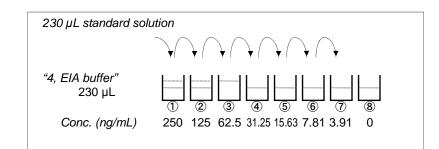

5) Dilution of test sample

Test samples should be diluted with "4, EIA buffer" as necessary. Serum samples need to be diluted with "4, EIA buffer" more than twice. If the concentration of human Napsin A in samples may not be estimated in advance, the pre-assay with several different dilutions will be recommended to determine the proper dilution of samples.

3. Measurement procedure

All reagents shall be brought to room temperature approximately 30 minutes before use. Then mix it gently and completely before use. Make sure of no change in quality of the reagents. Standard curve shall be prepared simultaneously with the measurement of test samples.

	Test Sample	Standard	Test Sample Blank	Reagent Blank
Reagents	Test sample 100 µL	Diluted standard (Tube 1-7) 100 μL	EIA buffer (Tube-8) 100 μL	EIA buffer 100 μL
Incubation overnight at 4°C with plate lid				
4 times (wash buffer more than 350 μL)*				
Labeled Antibody	100 µL	100 µL	100 µL	-
Incubation for 30 minutes at 4°C with plate lid				
5 times (wash buffer more than 350 μL)*				
Chromogen	100 µL	100 µL	100 µL	100 µL
Incubation for 30 minutes at room temperature (shielded)				
Stop solution	100 µL	100 µL	100 µL	100 µL
Read the plate at 450nm against a Reagent Blank within 30 minutes after addition of Stop solution.				

- 1) Determine wells for reagent blank. Put 100 μL each of "4, EIA buffer" into the wells.
- 2) Determine wells for test sample blank, test sample and diluted standard. Then, put 100 µL each of test sample blank (tube-8), test sample and dilutions of standard (tube-1-7) into the appropriate wells.
- Incubate the precoated plate overnight at 4°C after covering it with plate lid. 3)
- Wash the plate with the prepared wash buffer and remove all liquid.* 4)
- Pipette 100 µL of labeled antibody solution into the wells of test samples, 5) diluted standard and test sample blank.
- Incubate the precoated plate for 30 minutes at 4°C after covering it with plate 6) lid.
- Wash the plate with the prepared wash buffer and remove all liquid.* 7)
- Take the required quantity of "6, Chromogen" into a disposable test tube. Then, 8) pipette 100 µL from the test tube into the wells. Please do not return the rest of the test tube to "6, Chromogen" bottle to avoid contamination.
- Incubate the precoated plate for 30 minutes at room temperature in the dark. 9) The liquid will turn blue by addition of "6, Chromogen".
- Pipette 100 µL of "7, Stop solution" into the wells. Mix the liquid by tapping the 10) side of precoated plate. The liquid will turn yellow by addition of "7, Stop solution".
- Remove any dirt or drop of water on the bottom of the precoated plate and 11) confirm there is no bubble on the surface of the liquid. Then, run the plate reader and conduct measurement at 450 nm against a reagent blank. The measurement shall be done within 30 minutes after addition of "7, Stop solution".

SPECIAL ATTENTION

Specify the following concentration of each tube

c following co		
Tube-1	250 ng/mL	
Tube-2	125 ng/mL	
Tube-3	62.5 ng/mL	
Tube-4	31.25 ng/mL	
Tube-5	15.63 ng/mL	
Tube-6	7.81 ng/mL	
Tube-7	3.91 ng/mL	
Tube-8	0 ng/mL	(Test Sample Blank)

Put 230 µL of Standard solution into tube-1 and mix it gently. Then, put 230 µL of tube-1 mixture into tube-2. Dilute two times standard solution in series to set up 7 points of diluted standard between 250 ng/mL and 3.91 ng/mL. Tube-8 is the test sample blank as 0 ng/mL.


- Test samples should be measured soon after collection. For the storage ot test samples, store them frozen and do not repeat freeze/thaw cycles. Thaw the test samples at a low temperature and mix them completely before measurement.
- Test samples should be diluted with "4, EIA buffer", as the need arises. 2)
- Duplicate measurement of test samples and standard is recommended. 3)
- 4) Use test samples in neutral pH range. The contaminations of organic solvent may affect the measurement.
- 5) Use only wash buffer contained in this kit for washing the precoated plate. Insufficient washing may lead to the failure in measurement.
- 6) Remove the wash buffer completely by tapping the precoated plate on paper towel. Do not wipe wells with paper towel.
- "6, Chromogen" should be stored in the dark due to its sensitivity against light. 7) "6, Chromogen" should be avoided contact with metals.
- Measurement should be done within 30 minutes after addition of "7, Stop 8) solution".



CALCULATION OF TEST RESULT

Subtract the absorbance of test sample blank from all data, including standards and unknown samples before plotting. Plot the subtracted absorbance of the standards against the standard concentration on log-log graph paper. Draw the best smooth curve through these points to construct the standard curve. Read the concentration for unknown samples from the standard curve.

Example of standard curve

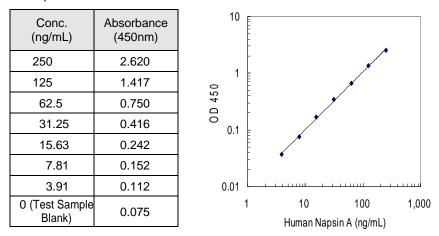

* The typical standard curve is shown above. This curve can not be used to derive test results. Please run a standard curve for each assay.

PERFORMANCE CHARACTERISTICS

1. Added Recovery Assay

Specimen	Theoretical Value (ng/mL)	Measurement Value (ng/mL)	%
10%FCS added RPMI-1640 (x2)	250.00	250.70	100.3
	125.00	114.92	91.9
	62.50	57.36	91.8
	31.25	26.79	85.7
	141.11	98.02	69.5
Human Serum	78.61	54.61	69.5
(x2)	47.36	34.09	72.0
	31.73	23.92	75.4
	131.52	99.78	75.9
Human Serum	69.02	49.89	72.3
(x4)	37.77	26.75	70.8
	22.15	16.11	72.7
	134.79	93.61	69.4
Human Urine	72.29	49.94	69.1
(x2)	41.04	28.86	70.3
	25.42	18.92	74.4
	129.93	100.17	77.1
Human Urine	67.43	51.32	76.1
(x4)	36.18	27.15	75.0
	20.55	15.35	74.7

2. Intra - Assay

Measurement Value (ng/mL)	SD value	CV value (%)	n
129.01	7.24	5.6	26
23.10	1.43	6.2	26
7.86	0.77	9.8	26

PRECAUTION FOR INTENDED USE AND/OR HANDLING

- 1. All reagents should be stored at 2 8°C. All reagents shall be brought to room temperature approximately 30 minutes before use.
- 2. "3, Standard" is lyophilized products. Be careful to open this vial.
- 3. "7, Stop solution" is a strong acid substance. Therefore, be careful not to have your skin and clothes contact "7, Stop solution" and pay attention to the disposal of "7, Stop solution".
- 4. Dispose used materials after rinsing them with large quantity of water.
- 5. Precipitation may occur in "2, Labeled antibody Conc.", however, there is no problem in the performance.
- 6. Wash hands after handling reagents.
- 7. Do not mix the reagents with the reagents from a different lot or kit.
- 8. Do not use expired reagents.
- 9. This kit is for research purpose only. Do not use for clinical diagnosis.

STORAGE AND THE TERM OF VALIDITY

Storage Condition : 2 - 8°C

The expiry date is specified on outer box.

REFERENCE

- Hirano T, Auer G, Maeda M, Hagiwara Y, Okada S, Ohira T, Okuzawa K, Fujioka K, Franzén B, Hibi N, Seito T, Ebihara Y, Kato H. Human tissue distribution of TA02, which is homologous with a new type of aspartic proteinase, napsin A. Jpn J Cancer Res. 2000 Oct;91(10):1015-21.
- Hirano T, Gong Y, Yoshida K, Kato Y, Yashima K, Maeda M, Nakagawa A, Fujioka K, Ohira T, Ikeda N, Ebihara Y, Auer G, Kato H. Usefulness of TA02 (napsin A) to distinguish primary lung adenocarcinoma from metastatic lung adenocarcinoma. Lung Cancer. 2003 Aug;41(2):155-62.
- Dejmek A, Naucler P, Smedjeback A, Kato H, Maeda M, Yashima K, Maeda J, Hirano T. Napsin A (TA02) is a useful alternative to thyroid transcription factor-1 (TTF-1) for the identification of pulmonary adenocarcinoma cells in pleural effusions. Diagn Cytopathol. 2007 Aug;35(8):493-7.

Version 2.

September 2016 *

Made in Japan.

3. Inter - Assay

Measurement Value (ng/mL)	SD value	CV value (%)	n
122.19	7.36	6.0	3
23.78	1.08	4.5	3
8.37	1.41	16.9	3

4. Sensitivity

2.92 ng/mL

The sensitivity for this kit was determined using the guidelines under the National Committee for Clinical Laboratory Standards (NCCLS) Evaluation Protocols. (National Committee for Clinical Laboratory Standards Evaluation Protocols, SC1, (1989) Villanova, PA: NCCLS.)

Immuno-Biological Laboratories Co., Ltd.